Angelcare®

FOR YOUR PEACE OF MIND

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TABLE OF CONTENTS

Warnings	4
Important safety instructions	4
Setting up your monitor	8
4 Easy Steps to Operate Your Angelcare® Monitor	9
Warnings	9
Instructions	9
Testing Your Angelcare® Monitor	10
Movement Sensitivity Adjustment	10
Troubleshooting	10
Warranty	11

E
N

Owner's manuel

KEEP THIS INSTRUCTION
MANUAL FOR FUTURE REFERENCE.

PLEASE READ THE WARNINGS AND
INSTRUCTIONS CAREFULLY BEFORE
ASSEMBLING & USING THE PRODUCT.

FAILURE TO FOLLOW THE WARNINGS
AND THE INSTRUCTIONS COULD RESULT
IN SERIOUS INJURY OR DEATH
OR DAMAGES TO THE PRODUCT.

⚠ WARNING

IMPORTANT SAFETY INSTRUCTIONS

THIS PRODUCT CANNOT REPLACE RESPONSIBLE ADULT SUPERVISION. This product is designed to aid in the monitoring of your child. When using the monitor, you must still provide the supervision necessary for the continued safety of your child. Monitoring of premature babies, or those considered to be at risk, should always be carried out under the supervision of a health professional. This monitor is not intended to be used as a medical monitor.

STRANGULATION HAZARD. Children have STRANGLED in cords. Keep the sensor pad cord out of the reach of children. Do not put the Nursery Unit inside a crib or play yard. Because of the serious strangulation risk to infants and toddlers, parents and caregivers should never place the unit within 3 feet (1 meter) of a crib.

Route sensor pad cord only below crib mattress and out of reach of children.

ALWAYS be sure to test the monitor before first use and then on a periodic basis.

ALWAYS be sure to immediately check on baby whenever an alarm sounds.

Do NOT use monitor near water (such as bath tub, sink, ect.).

Keep monitor AWAY from heat sources (such as stoves, radiators, ect.).

ALWAYS place the Nursery Unit on a flat surface in an upright position and allow for proper ventilation. Do NOT place on sofas, cushions, beds, etc. which may block ventilation.

This monitor uses public airwaves to transmit signals. The unit may pick up signals or interference from other nursery monitors in your area. Your

monitor may even be picked up by other households. To protect your privacy, make sure your monitor is turned off when not in use.

This product is NOT a toy. Do NOT allow children to play with it.

This product contains small parts. Adult assembly required. Exercise care when unpacking and assembling the product and keep small parts away from children.

Use and store the monitor at a temperature between 50°F (10°C) and 104°F (40°C).

⚠ BATTERY WARNING:

Do NOT mix old and new batteries.

Do NOT mix alkaline, standard or rechargeable batteries.

Insert batteries with correct polarity.

Do NOT short-circuit supply terminals.

ONLY use the recommended batteries or equivalent (volts and size).

Remove batteries when product is stored for long periods of time or when batteries are exhausted.

Dispose exhausted batteries properly.

Keep all batteries away from children.

⚠ IMPORTANT SAFETY INSTRUCTIONS

1) Read these instructions All the safety instructions should be read before the product is operated.

2) Keep these instructions. The safety and operating instructions should be retained for future reference.

3) Heed all warnings – All warnings on the product and in the operating instructions should be adhered to.

4) Follow all instructions – All operating and use instructions should be followed.

5) Water and moisture – Do NOT use this apparatus near water. For example, near a bath tub, wash bowl, kitchen sink, laundry tub, in a wet environment or near a swimming pool and the like. Care should be taken so that objects do not fall onto units and that liquids are not spilled into either the Nursery Unit or sensor pad through its openings.

6) Cleaning – Disconnect all units before cleaning. Do NOT immerse in water any part of the Nursery Unit. Clean only with dry cloth. Do NOT spray with cleaning products or solvents. Wipe sensor pad with a cloth moistened very slightly with an antiseptic or mild detergent.

7) Ventilation – Do NOT block any ventilation openings. Allow for proper ventilation when units are in use. Do NOT cover the Nursery Unit with any object such as a blanket. Do NOT place it in a drawer or in any location which would muffle the sound or interfere with the normal flow of air.

8) Heat sources – Do NOT install near any heat sources such as radiators, heat registers, stoves, TVs or other apparatus (including amplifiers) that produce heat. Heat can damage the case or electrical parts.

9) Power cord protection – Protect the sensor pad cord from being walked on or pinched by furniture or other items, paying particular attention to the point where they exit from the product.

10) Accessories – Only use attachments/ accessories specified by the manufacturer. Do NOT place this product on an unstable cart, stand, tripod, bracket or table. The product may fall, causing serious injury to a child or adult and serious damage to the product. Use only with a cart, stand, tripod, bracket or table recommended by the manufacturer or sold with the product. Any mounting of the product should follow the manufacturer's instructions.

11) Unplug this appliance during lightning storms or when not used for long periods of time.

12) Servicing – Do NOT attempt to service this product yourself as opening or removing covers (except battery covers) may expose you to dangerous voltage or other hazards.

13) Damage required servicing- Refer all servicing to qualified service personnel. Servicing is required when the apparatus has been damaged in anyway, such as cord is damaged, liquid has been spilled or objects have fallen into the apparatus, the apparatus has been exposed to rain or moisture, does not operate normally, or has been dropped.

14) Replacement parts - When replacement parts are required, be sure the service technician has used replacement parts specified by the manufacturer or have the same characteristic as the original part. Unauthorized substitution may result in fire, electric shock or other hazards.

15) Safety check- Upon completion of any service or repair to the product, ask the service technician to perform safety checks to determine that the product is in proper operating condition.

16) Usage-This baby monitor is intended as an aid. It is NOT a substitute for responsible and proper adult supervision and should not be used as such. Use this appliance ONLY for its intended use as described in this manual.

17) Object and Liquid Entry- NEVER push objects of any kind into this product through casing slots as they may touch dangerous voltage points or short out parts that could result in a risk of fire or electric shock. NEVER spill liquid of any kind on the product.

18) STRANGULATION HAZARD. Children have STRANGLED in cords. Keep the sensor pad cord out of the reach of children. Because of the serious strangulation risk to infants and toddlers, parents and caregivers should NEVER place the Nursery Unit within 3 feet (1 meter) of a crib.

19) Vibration Source- When the monitor is in use, do NOT use a crib mobile or other crib accessories that may produce vibrations. Other sources of vibrations may include fans, washing machines, loud music, etc. These sources of vibration may interfere with the performance of your Angelcare® monitor and the alarm will not sound if it continues to detect any movement other than your baby's. Be sure all vibration sources are eliminated before using the monitor.

You may experience false alarms for several reasons. The most common reason is the unit has not been turned off after you have removed your baby from the crib. Another reason is the Sensor Pad has not been placed under the mattress properly.

Be sure the Sensor Pad is on a flat, firm surface with the printed side facing up. The sensitivity level may also need to be adjusted. Refer to the Movement Sensitivity Adjustment section p.10. Always be sure to immediately check on baby whenever an alarm sounds.

20) Sensor Pad – Your sensor pad can be used in any crib that is stationary and designed to safely hold a baby. If the crib has a spring or slat base, place a hard board between the crib base and the Sensor Pad measuring at least 13" x 13" (33 x 33 cm) and 1/4 inch (6 mm) thick. Only use your Angelcare® monitor in a stationary sleeping environment such as a crib where the Sensor Pad can rest on a completely flat, firm surface.

ALWAYS route the sensor pad cord below crib mattress ONLY and out of reach of children.

IMPORTANT SAFETY INSTRUCTIONS RELATED TO USE OF BATTERIES

1) Use ONLY the type and size of batteries specified in the user manual. There is a risk of explosion if batteries are replaced with an incorrect type.

2) Do NOT dispose batteries in fire. They may explode. Check with local regulations for possible special disposal instructions.

3) Exercise care in handling batteries in order not to short out the batteries with conducting materials

such as rings, bracelets and keys. Overcharging, short circuiting, reverse charging, mutilation, or incineration of batteries must be avoided to prevent one or more of the following occurrences: release of toxic materials, release of hydrogen and/or oxygen gas, rise in surface temperature.

4) Do NOT attempt to recharge batteries identified for use with this product that are not intended to be charged. The batteries may leak corrosive electrolyte or explode.

5) Do NOT attempt to rejuvenate the batteries identified for use with this product by heating them. Sudden release of battery electrolyte may occur, causing burns or irritation to eyes or skin.

6) If and when necessary, the batteries must be replaced with identical new ones or equivalent (volts and size). If a battery to be replaced is connected with other batteries in series, replace the other batteries at the same time as well. Do NOT mix old and new batteries.

7) Reverse polarity installation of the batteries in the end product must be avoided. Reverse insertion of batteries can cause charging, and that may result in leakage or explosion.

8) Remove batteries from this product if you store the product over 30 days because the batteries could leak and damage the product.

9) Discard 'dead' batteries as soon as possible since 'dead' batteries are more likely to leak in a product. If a battery has leaked or vented, it should be replaced immediately using protective gloves.

10) Do NOT store this product, or the batteries identified for use with this product, in high-temperature areas. Batteries that are stored in a freezer or refrigerator for the purpose of extending shelf life should be protected from condensation during storage and defrosting. Batteries should be stabilized at room temperature prior to use after cold storage.

WARNING

Modifications not authorized by the manufacturer may void users authority to operate this device.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:


- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) This device must accept any interference received, including interference that may cause undesired operation.

Setting up your monitor

POSITION OF SENSOR PAD UNDER MATTRESS

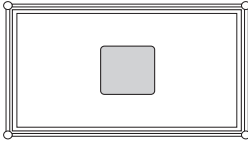MIDDLE OF THE CRIB

- Connect the Sensor Pad cord to the Nursery Unit.
- Place the Sensor Pad in the middle of the crib under the mattress, as illustrated besides. Make sure the printed side faces up. Do NOT place any bedding between the mattress and the sensor pad.
- Secure the sensor pad cord to crib using a tie wrap. **Run the cord on the floor to ensure it is out of the baby's reach.**
- The Sensor Pad must rest on a completely flat and firm surface. If the crib has a spring or slat base, place a hard board between the crib base and the Sensor Pad measuring at least 13" x 13" (33x33 cm) and 1/4 inch (6 mm) thick. Make sure the sensor pad's cord passes underneath the mattress as shown above.
- The sensor pad works with all baby bed sizes and all mattresses, except memory foam mattresses. You will have to test your sensor pad on STEP 1.
- Once you are done using the Sensor Pad, please remove it from the baby's crib and store it away.
- The Angelcare Sensor Pad feature is not designed for use with adults, children or infants older than 24 months.

WARNING

STRANGULATION HAZARD - Children have STRANGLED in cords. Keep the sensor pad cord out of the reach of children. Do not put the Nursery Unit inside a crib or play yard. Because of the serious strangulation risk to infants and toddlers, parents and caregivers should never place the nursery unit within 3 feet (1 meter) of a crib.

A Correct

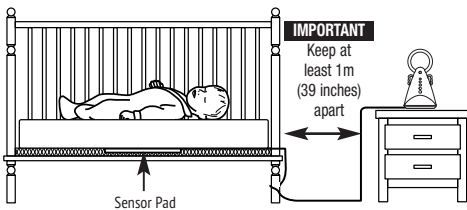
Figure 1

B Incorrect



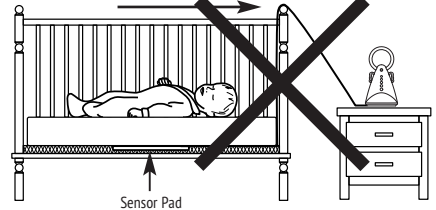
Figure 2

Figure 3

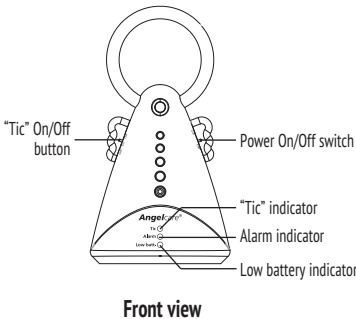
Front view

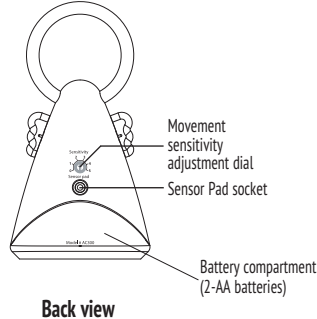Back view

Figure 4

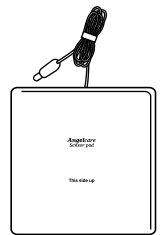Sensor Pad

4 EASY STEPS TO OPERATE YOUR ANGELCARE® MONITOR

STEP 1. SENSOR PAD

- Position Sensor Pad as illustrated (Figures 1 and 2). The printed side of Sensor Pad **MUST** face up. **DO NOT** tuck thick bedding under the edges of the mattress.
- The Sensor Pad must rest on a completely flat and firm surface. If the crib has a spring or slat base, place a hard board between the crib base and the Sensor Pad measuring at least 13" x 13" (33x33 cm) and 1/4 inch (6 mm) thick.
- The Sensor Pad works with all baby bed sizes and all mattresses, except memory foam mattresses or mattresses with a hollow frame. You will need to test your Sensor Pad after the set-up is completed.

STEP 2 MONITOR

- Insert 2 AA batteries into Monitor.
- Connect the Sensor Pad cord to the monitor (Figure 3).
- On the side of the Monitor, select 1 of the 2 functions of how you want to monitor your baby:
 "Tic" On: To see green movement indicator light flash **AND** hear a "Tic" sound with every movement baby makes.
 "Tic" Off: Same as above **BUT** without the "Tic" sound.

NOTE: The sensitivity adjustment has been preset. **DO NOT ADJUST UNLESS NECESSARY.** If poor capture of movement (red flashing light) or false alarm occurs, consult the Movement Sensitivity Adjustment section.

STEP 3 TO USE MONITOR

- Place baby in crib and switch Monitor to the ON position. When removing baby from crib, switch Monitor to the OFF position. If you forget to switch OFF, a "pre-alarm Tic" will sound after 15 seconds as a reminder.

STEP 4 BATTERY MODE

- The Angelcare® Movement Monitor is battery operated. Batteries will last approximately 3 to 4 months. When monitor produces a rapid clicking sound, change batteries promptly.

NOTE:

- Test your Angelcare® Monitor **BEFORE** first-time use and then on a periodic basis.
- When monitor is in use, **DO NOT** use crib mobile or other crib accessories that produce vibrations. Alarm will not sound as monitor is still detecting movement.
- The Angelcare® Monitor is perfectly safe. The Sensor Pad and the cord are entirely passive. They carry no electrical current and do not radiate any form of energy.
- The Sensor Pad works with all baby bed sizes and all mattresses, except memory foam mattresses or mattresses with a hollow frame. You will need to test your Sensor Pad after the set-up is completed.

Warnings

- This product is **NOT** a substitute for direct supervision of your baby. Check your baby's activities at regular intervals.
- The Sensor Pad may detect movement from many sources from **BOTH** inside and outside your baby's room, for example: fans, washing machines, loud music, etc. **Ensure that all vibration sources are eliminated before using the Angelcare® Monitor. Avoid contact with crib when monitor is in use. Refer to the Movement Sensitivity Adjustment and Troubleshooting sections.**
- **ALWAYS** place Monitor on flat surface in an upright position, away from baby's reach.
- This product should be serviced by qualified personnel. There are no user serviceable parts inside monitor. Refer all service questions to Angelcare.

Instructions

RETAIN INSTRUCTIONS: Keep user's manual for future reference.

WATER, MOISTURE AND HEAT DAMAGE: Keep all components away from water (sink, bath tub, swimming pool) and places that emit heat (stove, radiator).

VENTILATION: ALWAYS place Monitor on flat surface in an upright position so that air can flow freely around it.

CORD PROTECTION: Avoid damage to the Sensor Pad cord. The cord is long enough to route so that it is not likely to be walked on or pinched by items placed upon or against it.

CLEANING: Disconnect Monitor before cleaning. **DO NOT** immerse in water. Use a dry cloth to keep Monitor free from dust. **DO NOT** spray with cleaning products or solvents. Wipe Sensor Pad with a very slightly damp cloth with an antiseptic or mild detergent.

OBJECT AND LIQUID ENTRY: Care should be taken so that objects do not fall onto Monitor and that liquids are not spilled into either Monitor or Sensor Pad through its openings.

NON-USE PERIODS: When not using monitor for long periods, remove the batteries from Monitor.

Testing Your Angelcare® Monitor

NOTE: Test your Angelcare® Monitor **BEFORE** first-time use and then on a periodic basis.

The Angelcare® Monitor can be used outside of the crib, in a different sleeping area. Test monitor when moving Sensor Pad to a new location to ensure sensitivity is adjusted correctly to the new padding or mattress. Consult the Movement Sensitivity Adjustment section.

Step 1 Switch the power and “Tic” function to the ON positions. Gently place your hand on the mattress. The green light will continue to flash with your every movement. Remove hand from mattress. Due to no movement, a “pre-alarm Tic” sound will emit after 15 seconds. 5 seconds later, the alarm will sound and the red light will flash continuously.


NOTE: If alarm does not sound, the Sensor Pad may be detecting a continuous vibration from the floor, a strong draft or the crib is being touched. Avoid contact with crib when monitor is in use. Place crib near a solid supporting wall and away from strong air currents. The sensitivity of the Sensor Pad may need to be reduced, consult the Movement Sensitivity Adjustment section.

Step 2 To stop alarm, place your hand lightly on mattress for Sensor Pad to redetect movement. The red Alarm indicator light on the Monitor will continue to flash 1 minute (Figure 3).

Step 3 When Monitor produces a rapid “clicking” sound, replace batteries promptly.

Movement Sensitivity Adjustment

The Sensitivity Adjustment Dial has been preset to suit most situations. However in certain circumstances, the sensitivity may not be ideal. **DO NOT** increase sensitivity unless false alarms occurs and you have followed the advice in the Troubleshooting section for false alarms. Make a note of the pre-set position of the dial should you wish to return to it.

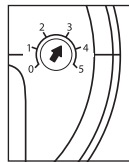The Alarm Indicator Light on the Monitor will flash red when Monitor barely detects movement. This could be a result of: baby being in a very deep sleep or baby having moved away from Sensor Pad to edge of crib. If alarm does **NOT** sound, no adjustment of sensitivity is required. If the Alarm Indicator Light flashes red most of the time with baby directly over the location of either Sensor Pad, a false alarm may occur. If this happens, adjust the sensitivity.

To adjust sensitivity, increase dial to the next setting. Retest monitor.

TROUBLESHOOTING

FALSE ALARMS

- Turn Monitor to the OFF position.
- Verify the connection between the Sensor Pad and Monitor. Unplug and reconnect the cord.
- Make sure there is no bedding between mattress and the Sensor Pad. The Sensor Pad must rest on a completely flat and firm surface. If the crib has a spring or slat base, place a hard board between the crib base and the Sensor Pad measuring at least 13” x 13” (33x33 cm) and 1/4 inch (6 mm) thick.
- Make sure Sensor Pad is positioned as described in the 4 Easy Steps section. Sensitivity may need to be adjusted, consult the Movement Sensitivity Adjustment section.

ALARM WILL NOT SOUND

- Avoid contact with the crib when monitor is ON.
- Place crib near a solid supporting wall to avoid Sensor Pad from detecting movement from outside crib.
- You may need to decrease the sensitivity, consult the Movement Sensitivity Adjustment section.

ALARM INDICATOR LIGHT IS FLASHING RED (MONITOR IS BARELY DETECTING MOVEMENT)

- Make sure there is no bedding between mattress and the Sensor Pad. The Sensor Pad must rest on a completely flat and firm surface. If the crib has a spring or slat base, place a hard board between the crib base and the Sensor Pad measuring at least 13” x 13” (33x33 cm) and 1/4 inch (6 mm) thick.
- Make sure Sensor Pad is positioned as described in the 4 Easy Steps section. Sensitivity may need to be adjusted, consult the Movement Sensitivity Adjustment section.

NO POWER

- Check / replace batteries.
- Make sure Monitor is ON.

LIMITED WARRANTY

Angelcare® warrants to the original owner of this product that the Angelcare® Movement Monitor (Model AC300) is free from defects in materials and workmanship for a period of one (1) year from the date of the original purchase. If this monitoring system fails to function properly when used as directed under normal conditions, within the one year period, Angelcare® will either repair or replace the product at our option free of charge, when it is returned complete and in its original packaging, prepaid postage, and insured to Angelcare.

The product must be accompanied by a proof of purchase, either a bill of sale, receipted invoice, or other proof that the monitoring system is within the warranty period. Angelcare® will bear the cost of repairing or replacing the product and shipping it back to you.

This warranty does not apply to a product that has been damaged as a result of improper maintenance, accidents, improper voltage supply or other misuse. The warranty is also void if the owner repairs or modifies the product in any way. Angelcare® is not liable for any incidental or consequential damages with regard to this product. The warranty also excludes any liability other than what is stated above. No other warranties are given.

LEGAL RIGHTS VARY FROM COUNTRY TO COUNTRY. SOME COUNTRIES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATIONS MAY NOT APPLY TO YOU.

WE ARE NOT IN ANY WAY RESPONSIBLE FOR ANY OUTCOME RESULTING FROM USE.

IMPORTANT! This Angelcare® Monitor is a personal care product. Please do not return this product to the store. If you have questions or need help, please call Angelcare.

REGISTER YOUR ANGELCARE PRODUCTS

Register your Angelcare monitor and other Angelcare products on our website: www.angelcare-monitor.com

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IMPORTANT! This Angelcare® monitor is a personal care product.
Please do not return this product to the store.
If you have questions or need help, please call Angelcare.

**Replacement Parts
Warranty Information (USA & Canada only):**

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Protect the environment!

This electronic product contains reusable material and must not be disposed in domestic waste. Please take it to a collecting point for electric and electronic equipment (ask your local authority for details).

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B *I* U | | Format | Font size | |

Why this monitor?

- Can be used with the sound monitor you already own
- You already own the Angelcare® Movement & Sound Monitor but have more than one baby to monitor
- You feel you don't need a sound monitor since you live in small accommodations or baby sleeps in the same room
- Can be used anywhere baby sleeps (camping, chalet and other places without electricity)

Features

- Under-the-mattress Sensor Pad detects every movement baby makes
- Alarm sounds after 20 seconds if absolutely no movement is detected
- Adjustable sensitivity
- Battery operated on 2 AA cells

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AC300_fra.jpg
no file selected

General infos

Stock*
Included shipping
each*
Weight*
Weight unit*
Length*
Width*
Height*
Dimension unit*

> Save

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The ANGELCARE trademark is a registered trade mark of the Angelcare Development Inc.